
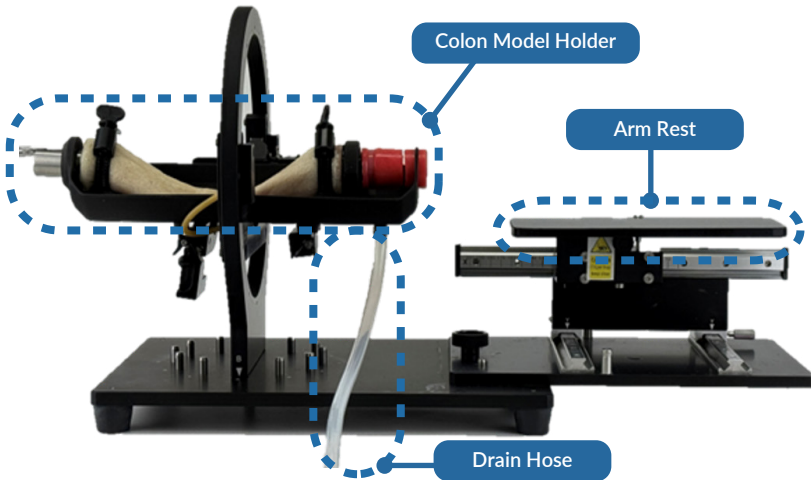


G-Master ESD Colon Tube Kit


PRODUCT INFORMATION

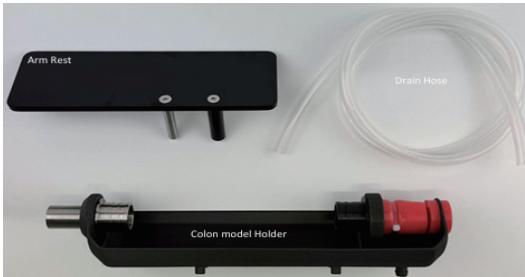

PRODUCT SPECIFICATIONS

Product Name	G-Master ESD Colon Tube Kit
Product Number	GM-C-TUBEKIT


DIRECTIONS FOR USE

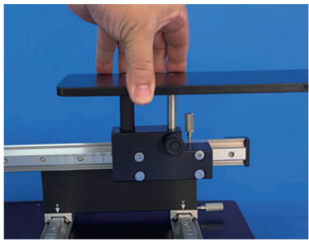

Detach the endoscope insertion section of the G-Master and install the hand rest.




Set the G-Master clip arms [L] and [R] to position ▼2.

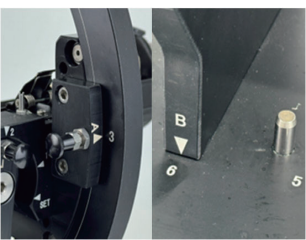

Set the G-Master [A] to [3] and [B] to [6].




Insert the two hooks on the underside of Colon Tube Fixation Part into the clips and fix them.

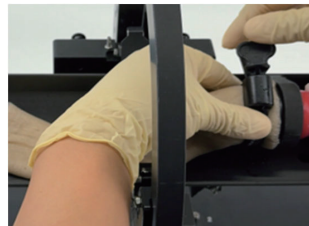

Set the Colon Model (sold separately) and fix it with the supplied band.



Set the Drain Hose to the hose port on the underside and direct the drainage into a bucket or similar container.
When performing irrigation, use normal saline or a saline solution with an equivalent salt concentration.


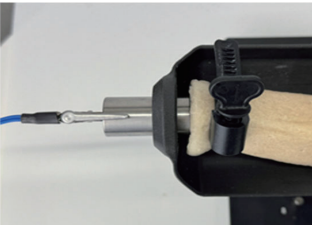


Insert the two hooks on the underside of Colon Tube Fixation Part into the clips and fix them.

Colon Model for ESD

PRODUCT INFORMATION

The **VTT Colon Model for ESD** is compatible with the existing G-Master using the **ESD Colon Tube Kit**, this model allows up to five simulation sessions and can be opened flat after incision for immediate visual confirmation of the results.

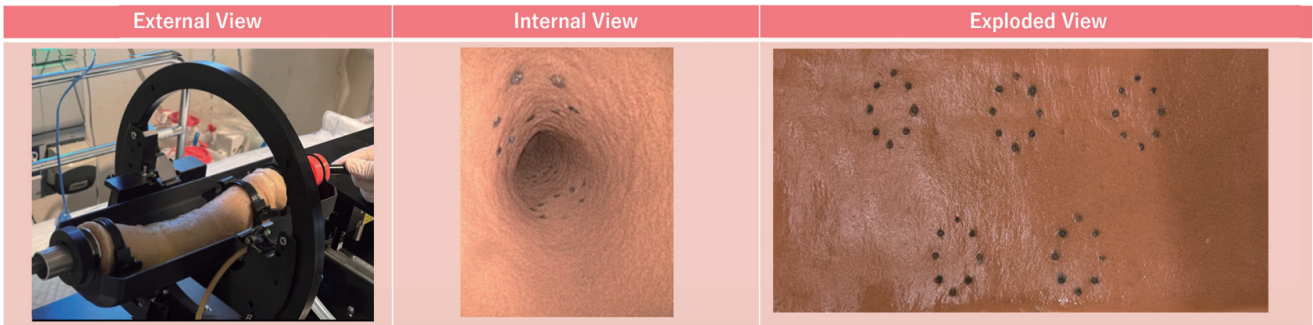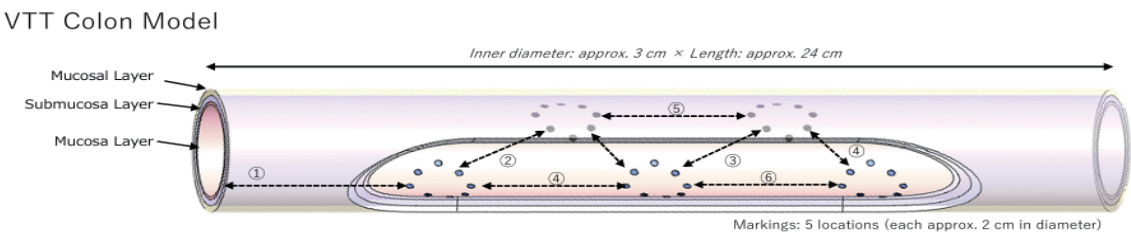

PRODUCT SPECIFICATIONS

Product Name	Colon Model for ESD
Product Number	VTT-ESD-TUBE01



ESD Colon Tube Kit

Designed for users who only need a colon model holder without the G-Master system.

PRODUCT SPECIFICATIONS

Product Name	ESD Colon Tube Kit
Product Number	GM-C-TUBEKIT-V2



Instruction Guide
&
Product Overview Video