


G-Master ESD Colon Tube Kit


User Manual

Product Number: GM-C-TUBEKIT

The G-Master main unit is sold separately.

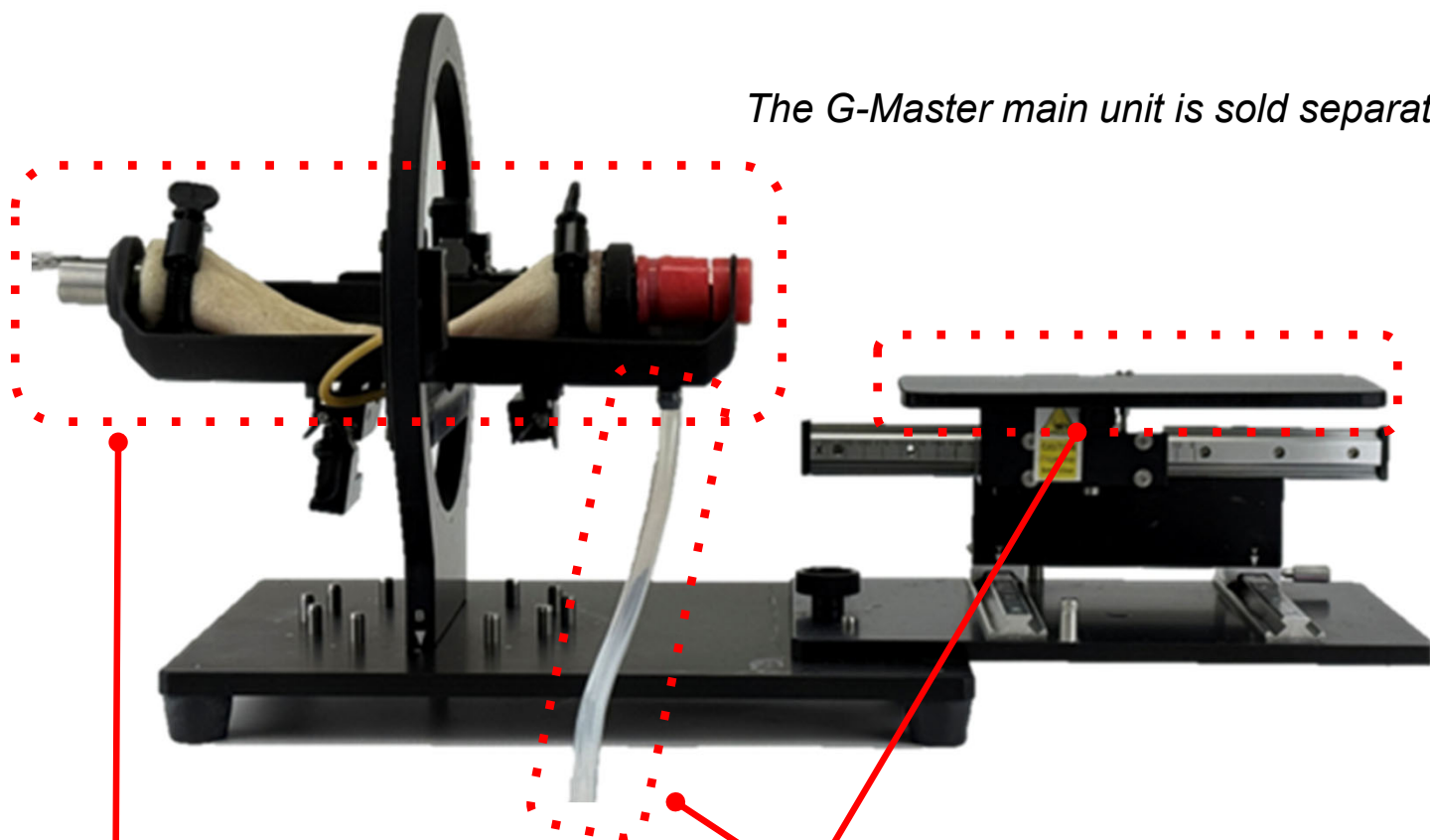

Product Components

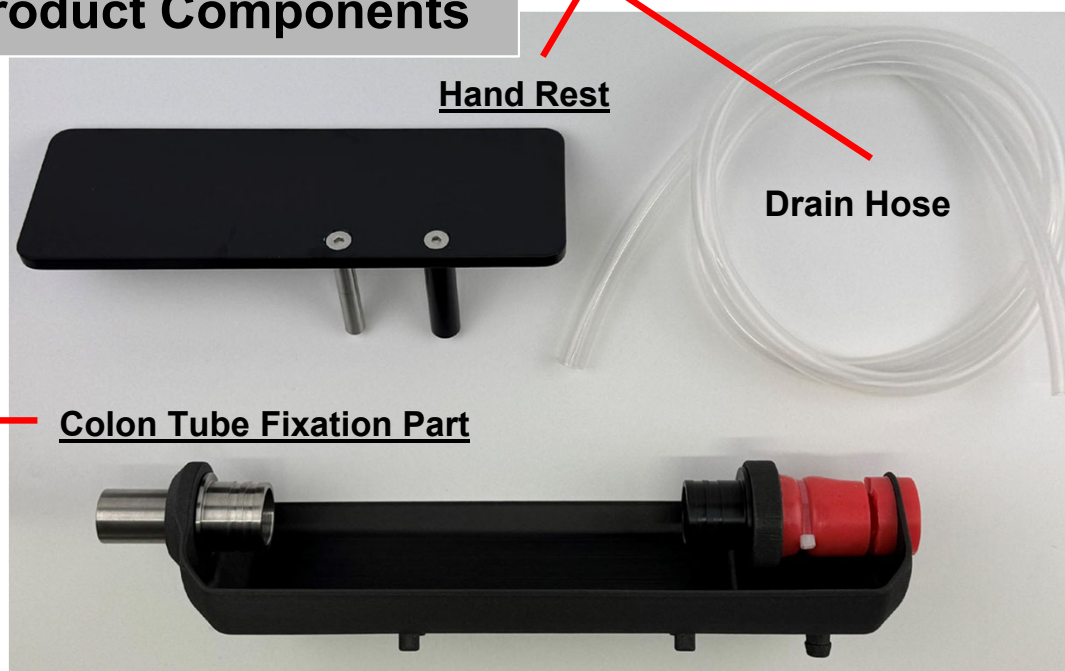

Hand Rest

Drain Hose

Colon Tube Fixation Part

KOTOBUKI Medical, Inc.

4-7-3 Chuo, Yashio, Saitama 3400816 Japan

<https://kotobukimedical.com/en/>

TEL +81 (0)48 951 5211

✉ info2@kotobukimedical.com






Setup And Directions


Scan the QR code to watch the setup video.



1. Detach the endoscope insertion section of the G-Master and install the hand rest. 2. Set the G-Master clip arms [L] and [R] to position ▼2. 3. Set the G-Master [A] to [3] and [B] to [6].

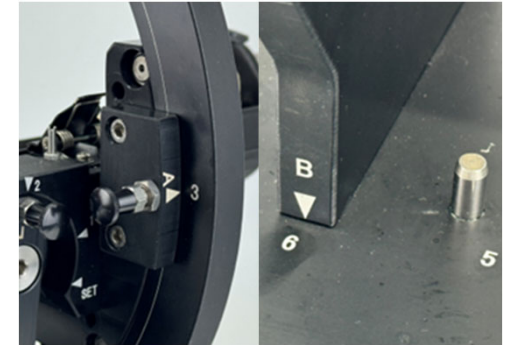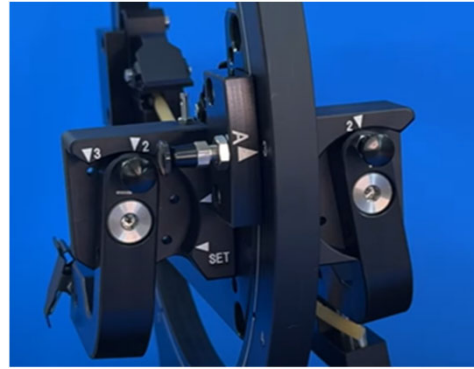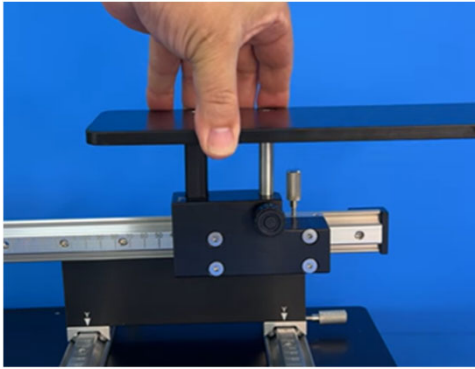

4. Insert the two hooks on the underside of Colon Tube Fixation Part into the clips and fix them. 5. Set the ESD Colon Tube Model (sold separately) and fix it with the supplied band.

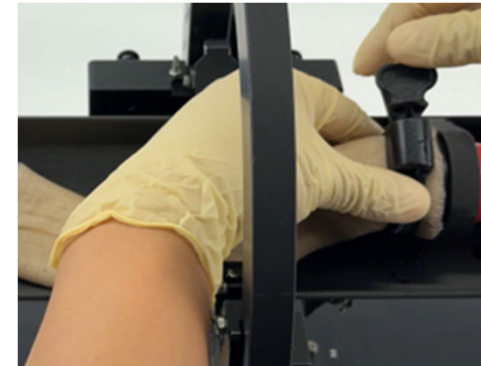

6. Set the Drain Hose to the hose port on the underside and direct the drainage into a bucket or similar container. 7. Connect the VTT Counter-Polar Cable to the metal tube at the end.

※When performing irrigation, use normal saline or a saline solution with an equivalent salt concentration.

