


# 11 posterior wall of the antrum

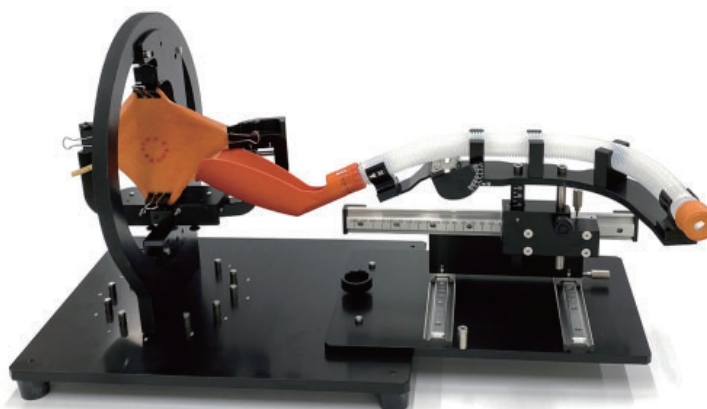

<https://kotobukimedical.com/en/esd-training/G-master/setting/>

\*For smooth setting, return the [X: Tube Left / Right] to the 0 position before setting.

\*Be careful not to pinch your fingers when adjusting the [Z: Tube Up / Down].

Mark	Name	Location Number											
<b>R</b>	Mucosa Holding Arm Right	SET		1 (straight)		2 (curve)		S : 1 (gassing) ⇔ 2 (degassing)					
<b>L</b>	Mucosa Holding Arm Left	SET		1 (straight)		2 (curve)		3 (degassing)					
<b>A</b>	Vertical Rotation position	1	2	3	4	5	6	7	8	9	10	11	12
<b>B</b>	Horizontal Rotation Position	1	2	3	4	5	6	7	8	9	10	11	12
<b>Y</b>	Tube Front / Rear		+20	+10	0	-10	-20	-30	-40	-50			
<b>X</b>	Tube Left / Right	0	10	20	30	40	50	60	70	80	90	100	
			110	120	130	140	150	160	170	180	190	200	210
<b>Z</b>	Tube Up / Down	0	10	20	25	30	35	40	45	50	60	70	80
<b>W</b>	Tube Angle	-20	-10	0	5	10	15	20	25	30	35	40	40
<b>M</b>	Spatula Insertion Position	OFF	1	2	3	4	5						
	Back Rubber Strap	OFF	ON										